

ToxGenie

Quick Installation & User's Manual


ToxGenie.com

Quick Installation & Windows Security Note

ToxGenie is a premium, GLP-compliant statistical software trusted by toxicology professionals. To provide this specialized tool at an accessible subscription rate rather than passing on the high costs of commercial publisher certificates, ToxGenie is distributed directly to our users.


As a result, Windows Defender SmartScreen may initially flag the installer.

SmartScreen Scenario 1:

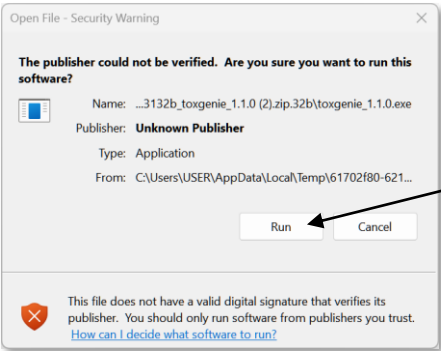


Click [Run]

SmartScreen Scenario 2:

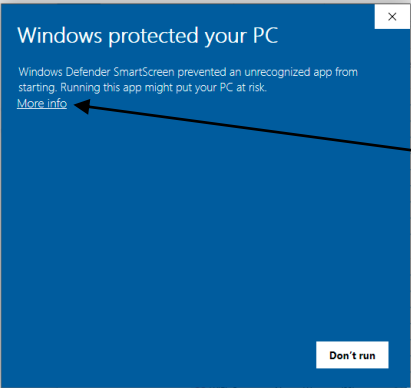

Step 1: On the blue screen, click [More info].

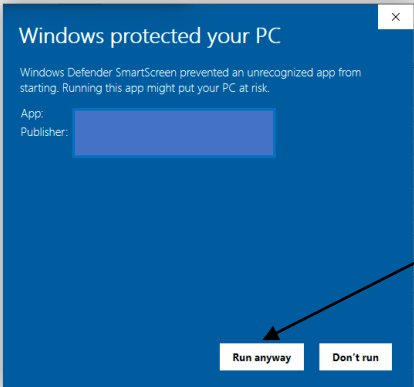

Step 2: Click [Run anyway] to proceed with the installation.

License Activation (Online & Offline)

Upon first launch, ToxGenie requires registration to activate your 30-day free trial or paid subscription. the installer.

Standard Activation (Online):

The screenshot shows the 'ToxGenie License Activation' window. At the top, it says 'Activate Your ToxGenie License'. There is a button labeled 'Get Free 30-Day License Key' (1). Below it is a text input field for 'Your Institutional Business Email:' (2). Underneath is a 'License Key:' section with a note 'Please copy and paste the entire License Key from your email.' and a corresponding text input field (3). At the bottom, there is an 'Activate License' button (4) and a 'Manual Offline' link.

Click **[Get Free 30-Day License Key]** to receive your key via email, or use your purchased key.

Enter your Institutional Business Email and the License Key.

Click **[Activate License]**.

Offline Mode (Institutional Firewall):

If your institution's strict firewall blocks the activation server, ToxGenie will automatically switch to **Offline Mode**.

The screenshot shows the 'ToxGenie License Activation' window in 'Offline Mode'. It states 'Online validation is blocked by your institution's firewall.' and provides instructions: '1. Email your Request Code below to toxgenie@toxgenie.com using your institutional email.' and '2. Copy and paste the License Key you receive from us into the field below.' There is a text input field for 'Your Request Code:' containing 'REQ-C8F63236' (1). Below it is a text input field for 'Enter your License Key:' (2). At the bottom, there is an 'Activate Offline' button (3) and a 'Cancel' button.

Copy the generated **Request Code**. and Email it to toxgenie@toxgenie.com.

Enter the manual offline key we provide to unlock the software.


Click **[Activate License]**.

⚠ Important: Institutional Email Required

Personal or commercial email accounts (e.g., @gmail.com, @yahoo.com) are **strictly prohibited**. You must use your official **institutional or business email address** (e.g., .edu, .ac, corporate domains) to receive a valid license key.

License Management & Premium Features

ToxGenie provides a transparent way to manage your subscription and access advanced GLP-compliant tools.

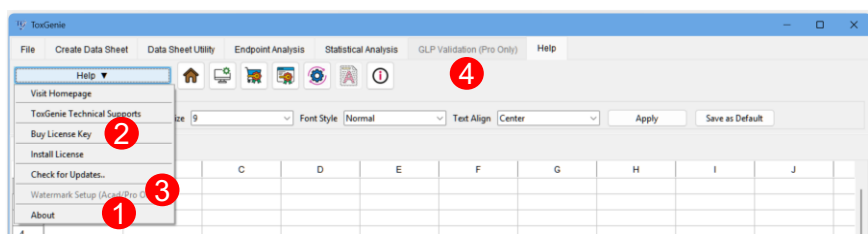

1 Checking License Validity:

To check your current subscription tier (Trial, Academic, Professional) and expiration date, navigate to the **[Help]** tab and click the **[About]** button.



2 Purchasing & Upgrading:

To purchase a new license or upgrade your current plan, click **[Buy License Key]** in the **Help** tab to visit our secure pricing page. Once you receive your key, click **[Install License]** to instantly activate it without losing your data.

Purchasing & Upgrading:


- 3 **Custom Watermarks (Acad & Pro):** Embed your official institutional watermarks and headers directly into the exported PDF and MS Word reports.
- 4 **GLP Validation Kit (Pro Only):** Unlock the **[GLP Validation]** tab to automatically generate out-of-the-box 21 CFR Part 11 compliance documentation, including complete IQ/OQ Protocols.

💡 **Note on Email Delivery:** Due to strict institutional security filters, the automated license key email may be routed to your **Spam or Junk folder**. Please check there if you do not receive it within a few minutes.

Quick Start: The "Zero-Manual" Workflow


ToxGenie is designed with a core philosophy: **Eliminate thick, complicated manuals.** If you understand toxicology, you already know how to use ToxGenie. Just follow the tabs from left to right.

Step 1: Prepare Data (Create Data Sheet Tab):

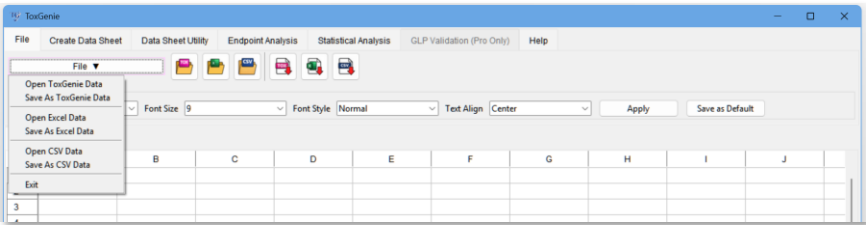

No need to struggle with data formatting. Simply go to the **[Create Data Sheet]** tab and click your study type (e.g., Quantal Data for OECD TG 203, Fish, Acute Test Data). A pre-formatted standard template will open instantly. (Alternatively, use the [File] tab to import your **Excel/CSV**

Step 2: Run Analysis (Endpoint / Statistics Tab)

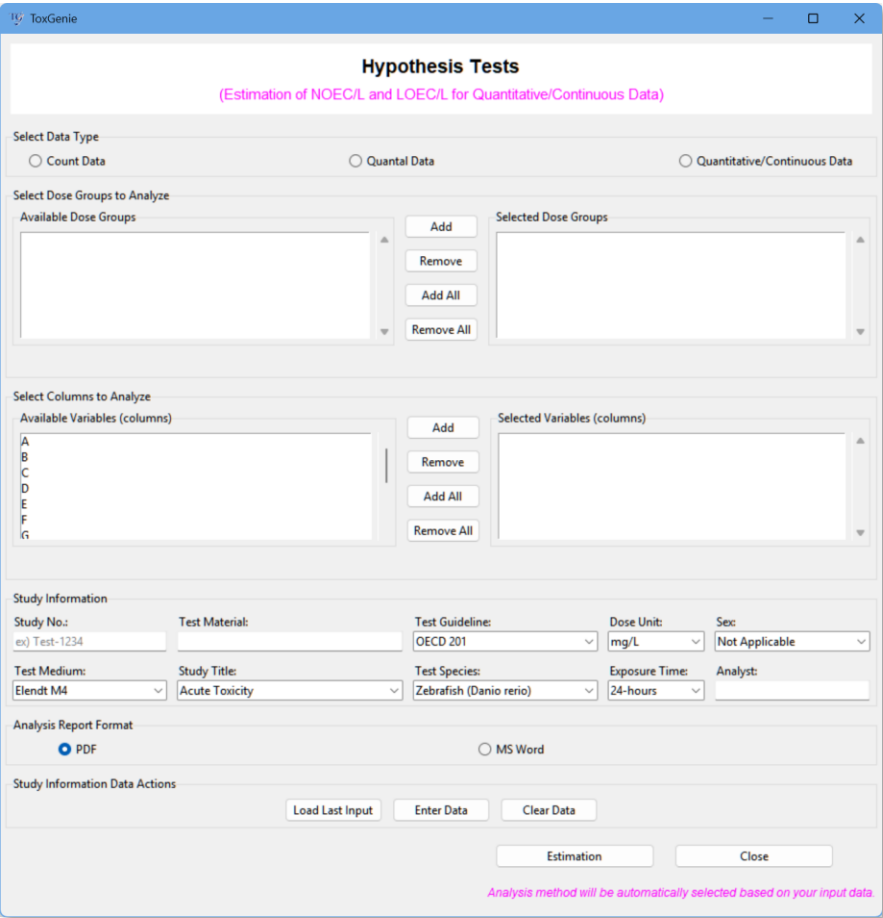

Quick Start: The "Zero-Manual" Workflow

Navigate to the **[Endpoint]** or **[Statistical Analysis]** tab and select a method (e.g., Hypothesis Tests). In the analysis window:

- Select your Data Type.
- Move target Doses and Variables to the 'Selected' boxes.
- Enter Study Information and choose the report format (PDF/Word).

Step 3: Generate Report

Click **[Estimation]**. Validated reports and graphs are generated instantly.

Critical Input Rules for Analysis

To ensure statistical integrity and avoid errors, ToxGenie enforces strict GLP validation rules:

- ✓ **Study No.:** Must be a combination of letters and numbers (e.g., Test-123). Entering numbers only is not allowed.
- ✓ **Quantal Data:** When selecting the Quantal Data type, you must specify the Initial Exposed Data Column from the dropdown list.
- ✓ **Dose Groups & Error Messages:** Depending on the specific statistical method you choose, non-numeric dose groups (e.g., "Positive Control") may or may not cause an error. **If an error or warning message box appears during estimation, please read the instructions carefully and follow them** (e.g., removing the non-numeric group from the Selected Dose Groups list if prompted).